



GUIDELINE OF BLZ CLOUD PLATFORM

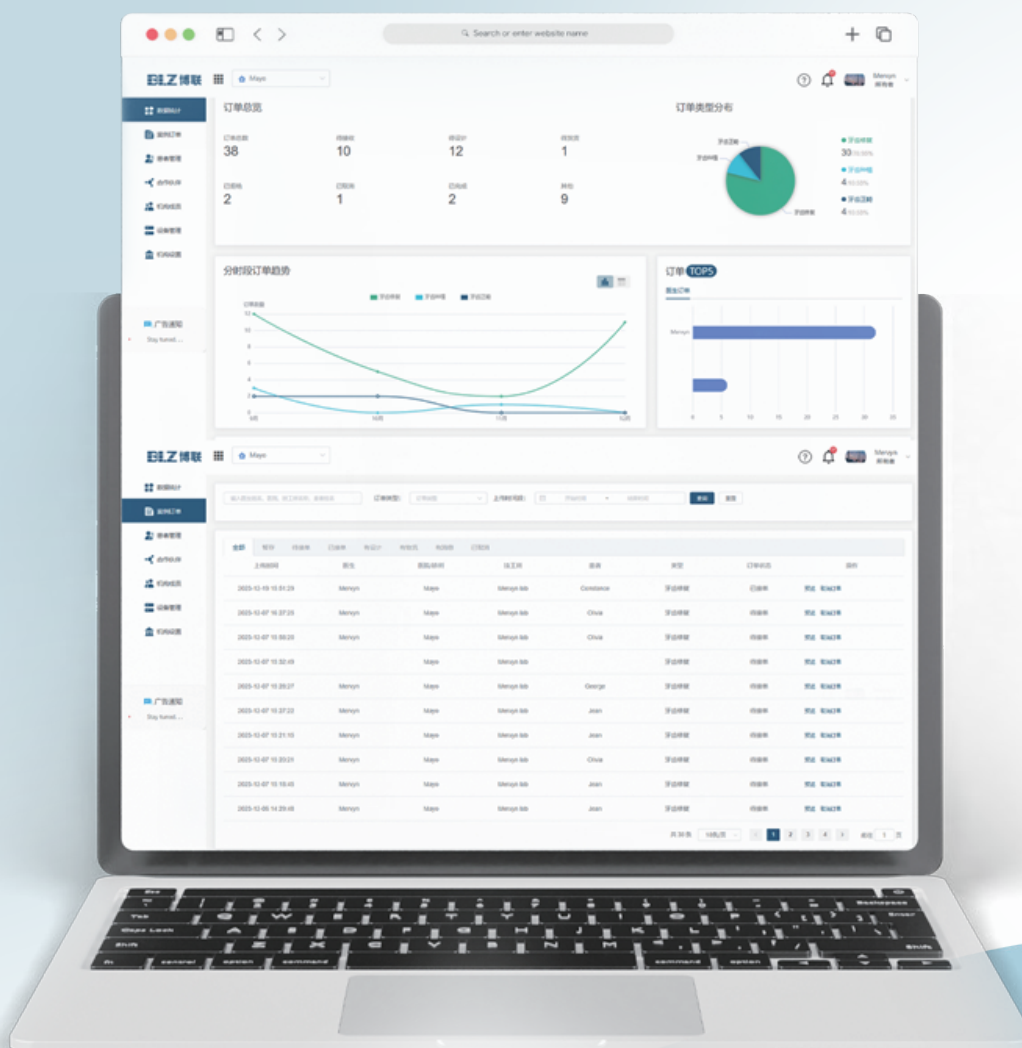

Content

- 1.Registration for a new user 1
- 2.Cooperation with a partner 2
- 3.Order management 9


User Flow for BLZ Cloud Platform

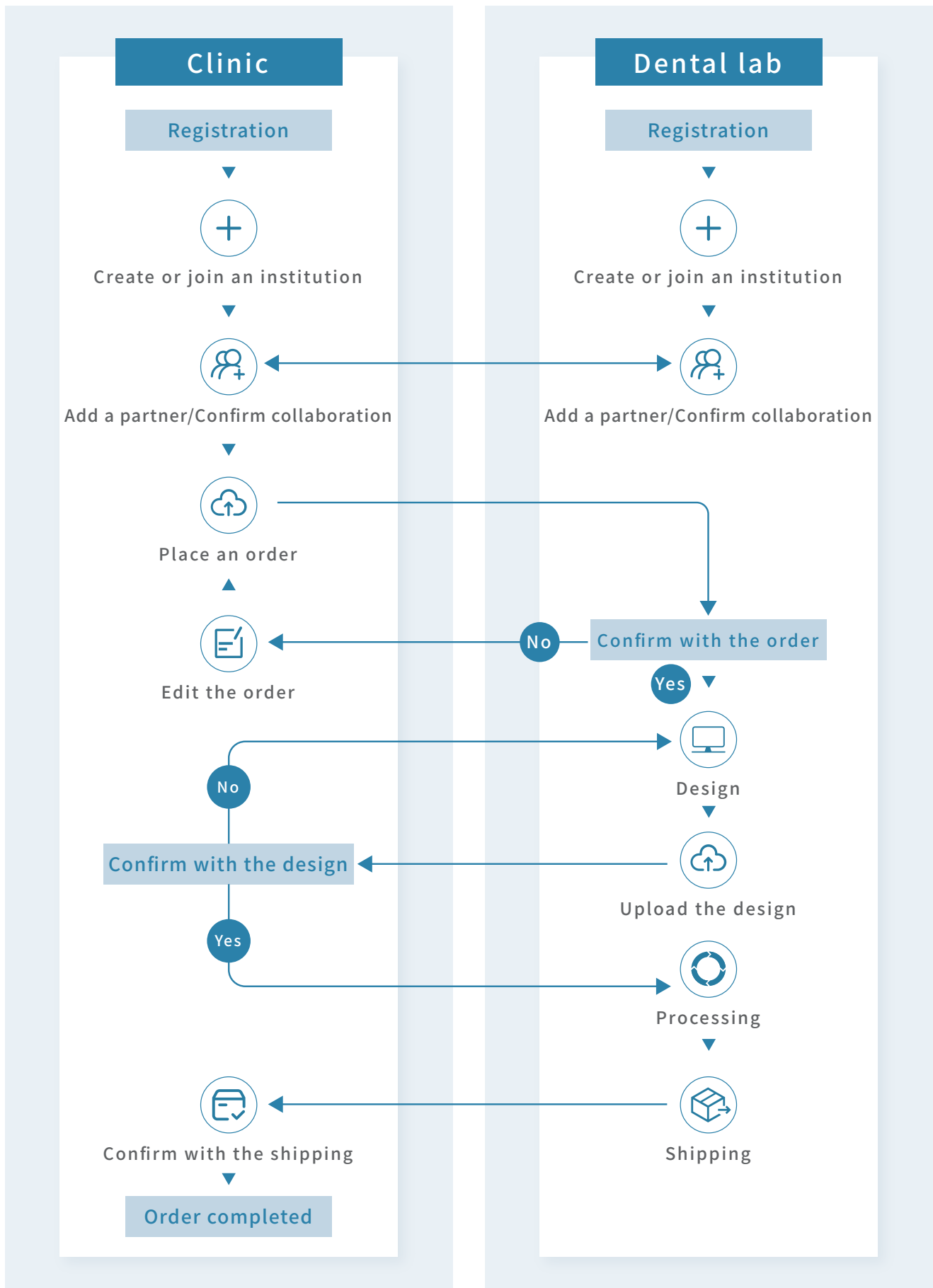

1.Registration for a new user

1.1 Registration and sign-in

- Enter the cloud site through software or copy the URL below to the browser
<https://www.blzdentalcloud.com/#/login>

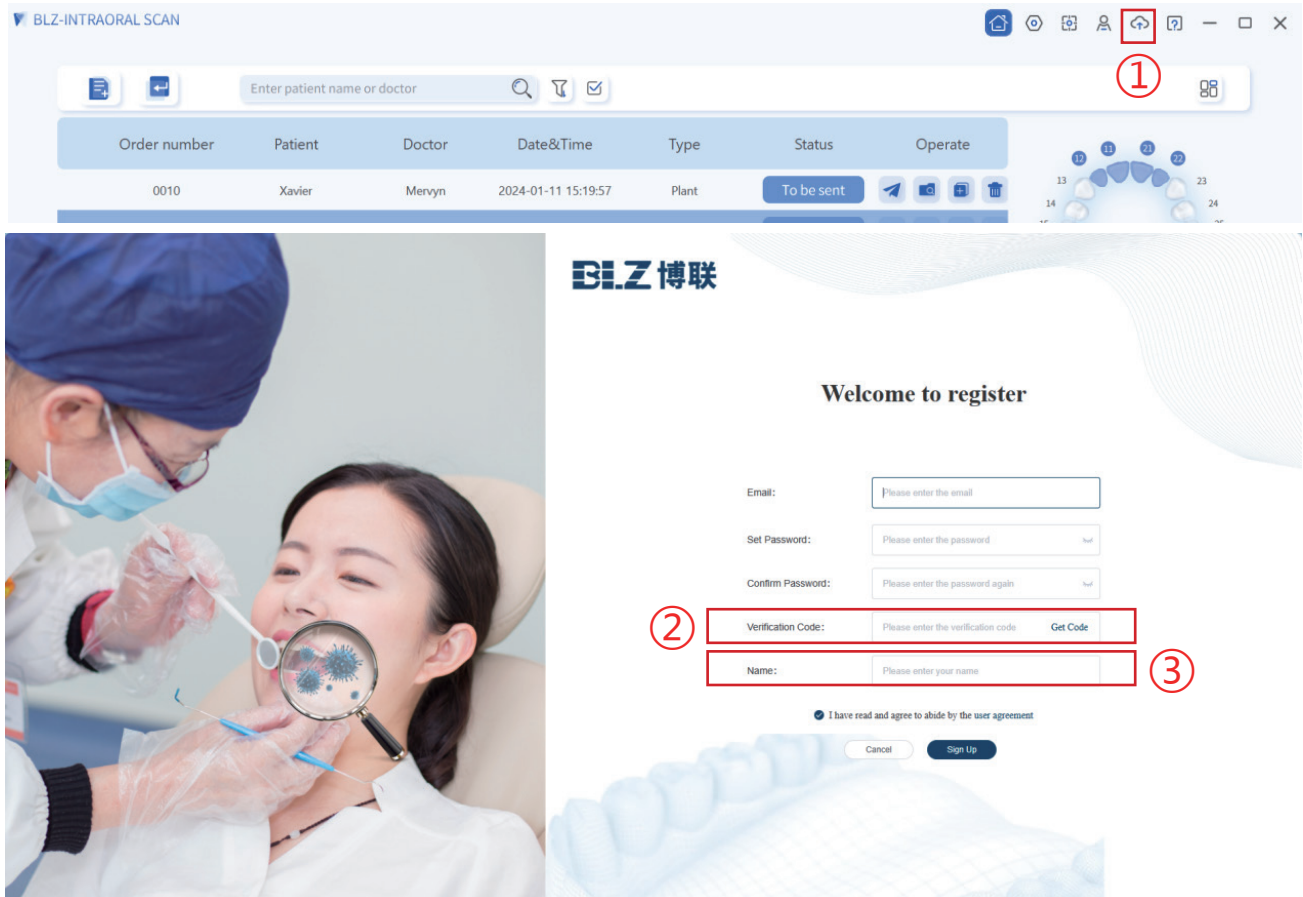


- ① Click to enter BLZ Cloud Platform. Enter registration interface by clicking “register”
- ② Verification code is sent to your registration email
- ③ Input name of your clinic or dental lab

- After registration, create a new institution or join an existing one to use the cloud normally.



- ① Click to join an existing clinic or dental lab
- ② Click to create your own clinic or dental lab

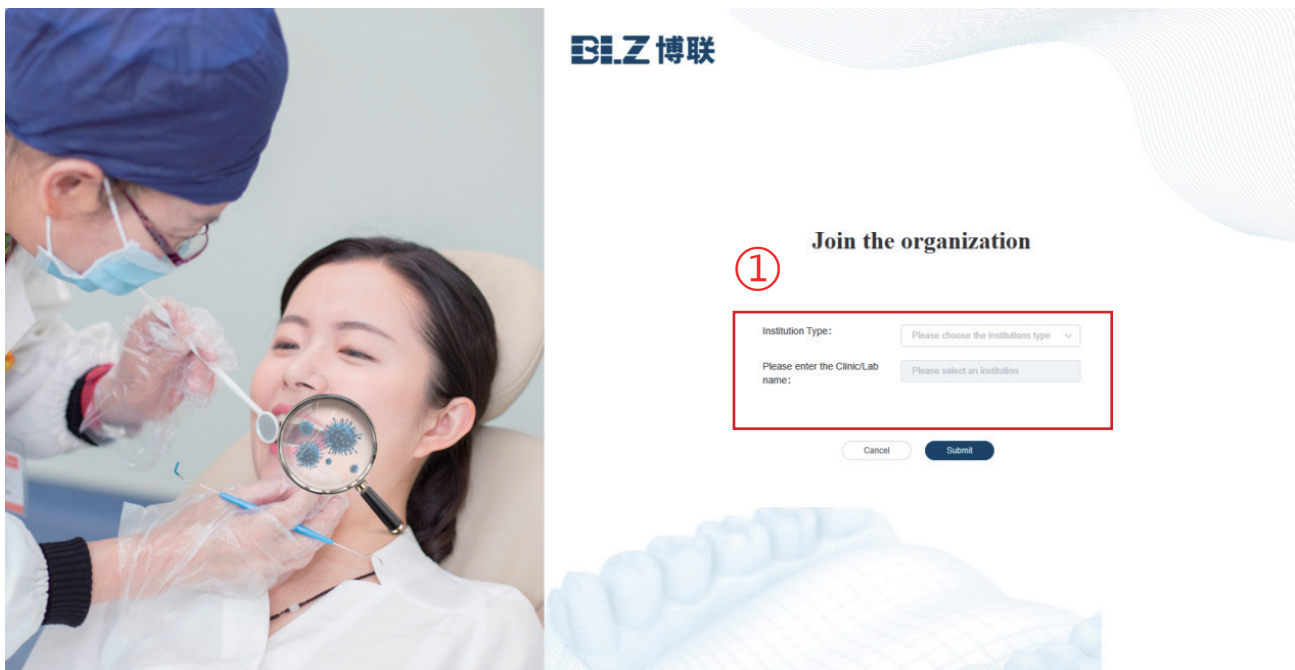

- ① Choose the institution type and input the name to join one




① Choose the institution type and input name to create one

1.2 New member and institution

1.2.1 New member

- Add a new member on institution side:

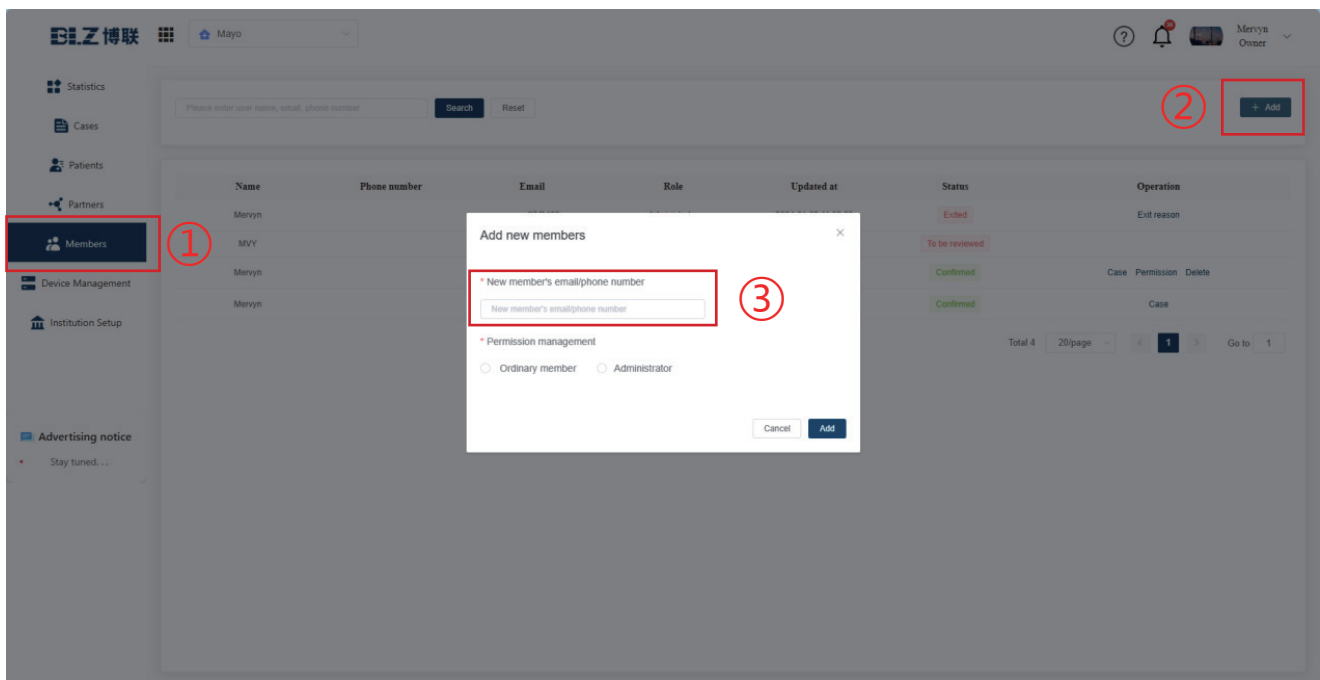


① Click “Members”

② Click “Add”


③ Input the member's email (register beforehand)

- The owner can manage members who've already joined.

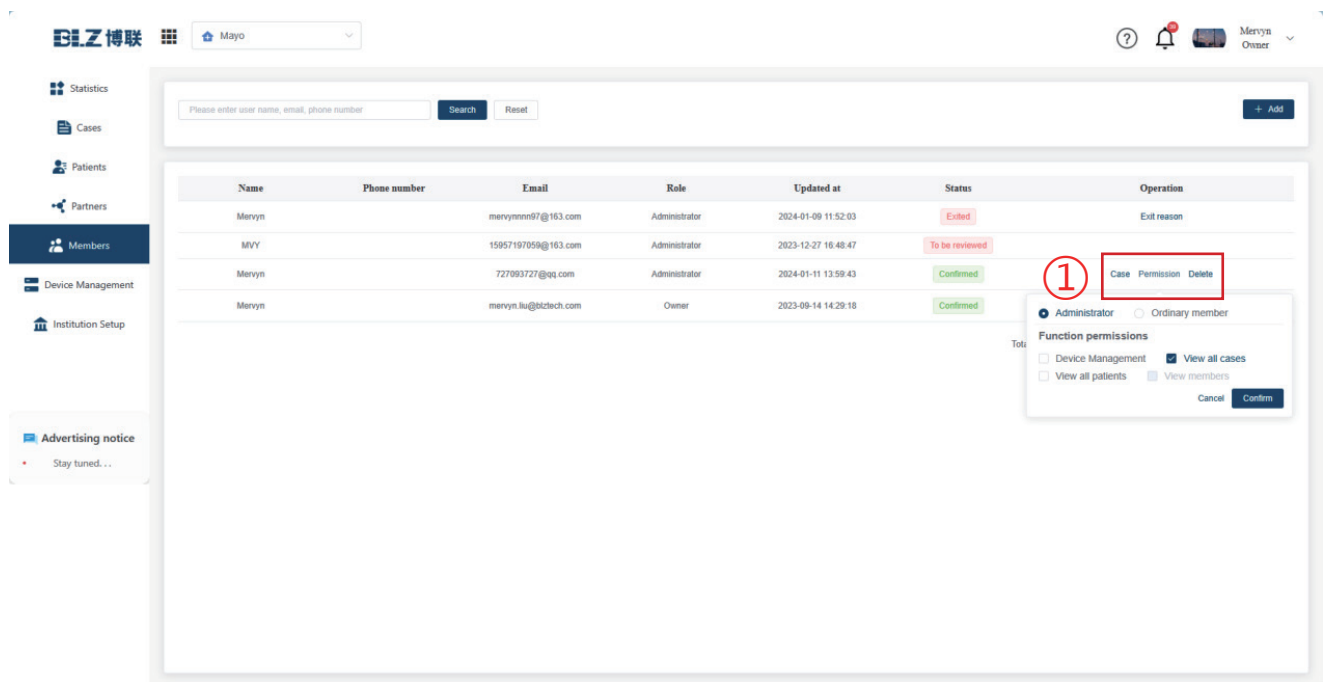


- ① Owner can set permissions for members

- New member logs in to agree the request.

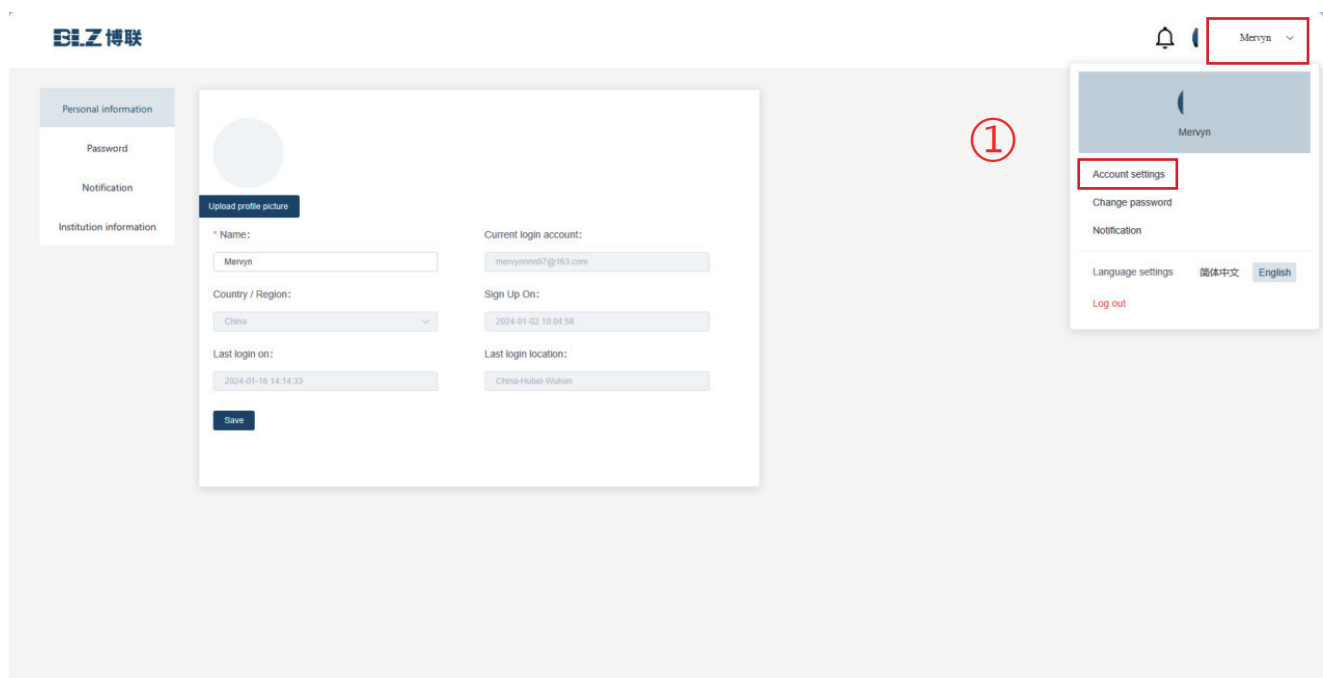

- ① Click "profile" - "Account settings"

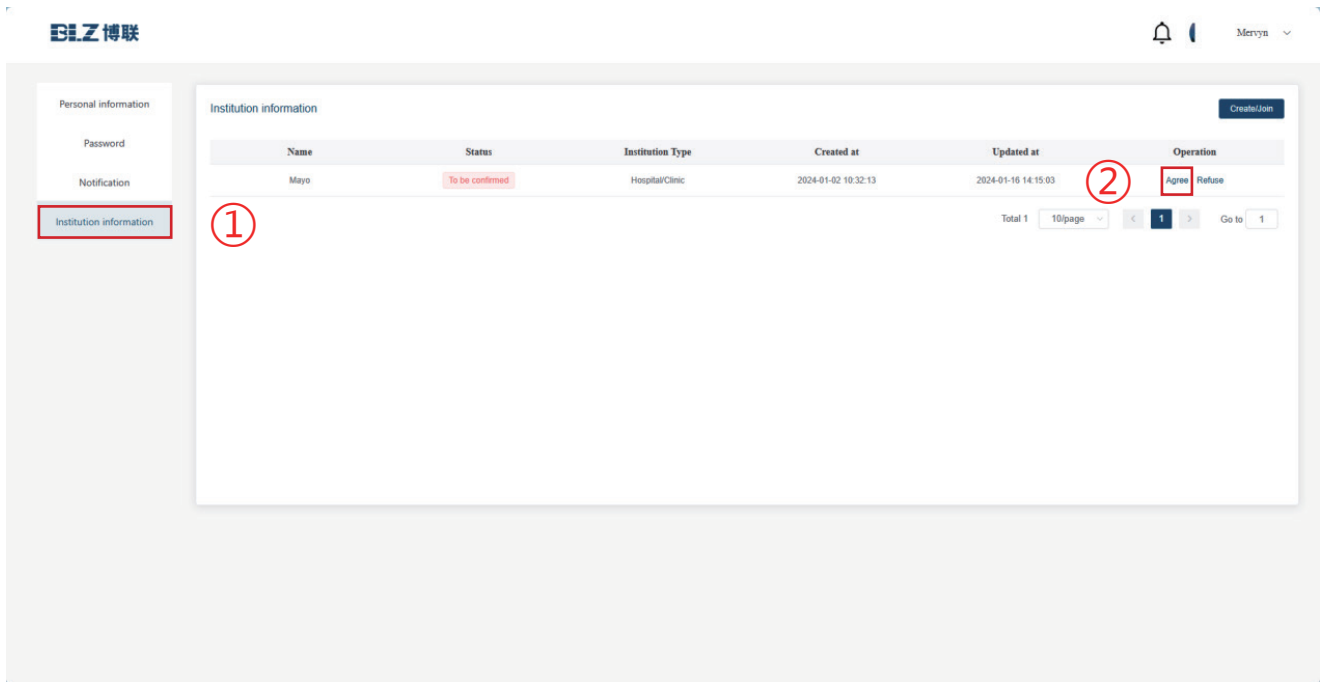


① Click “institution information”


② Choose “agree”

- New members can manage orders after obtaining relevant permissions.

1.2.2 New institution

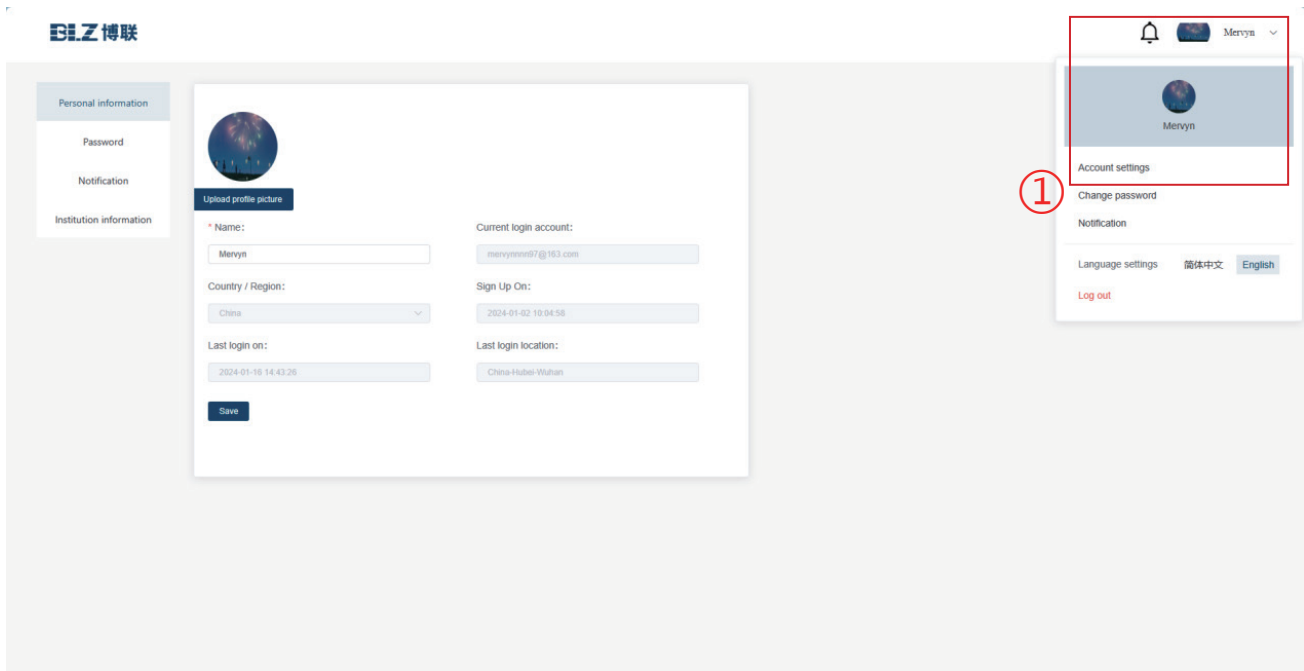


① Click “profile” - “Account settings” - “Institution information”

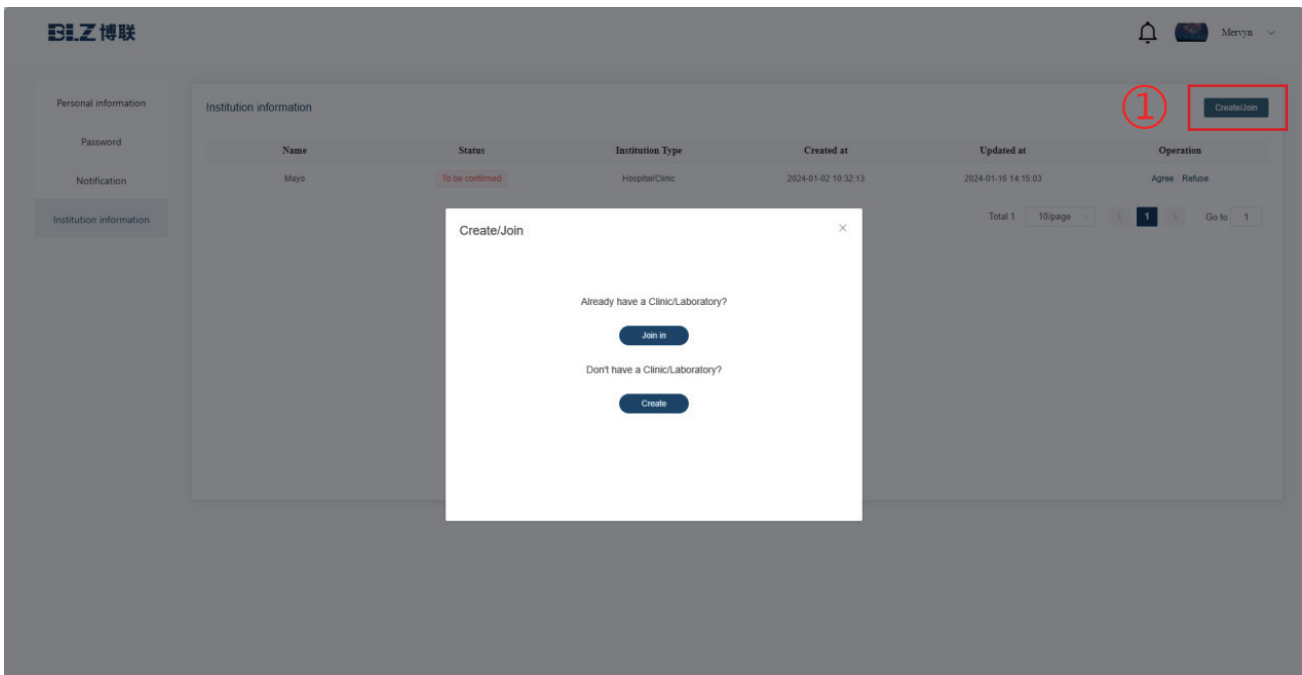


- ① Click “create” or “join” to create or join an institution

2.Cooperation with a partner

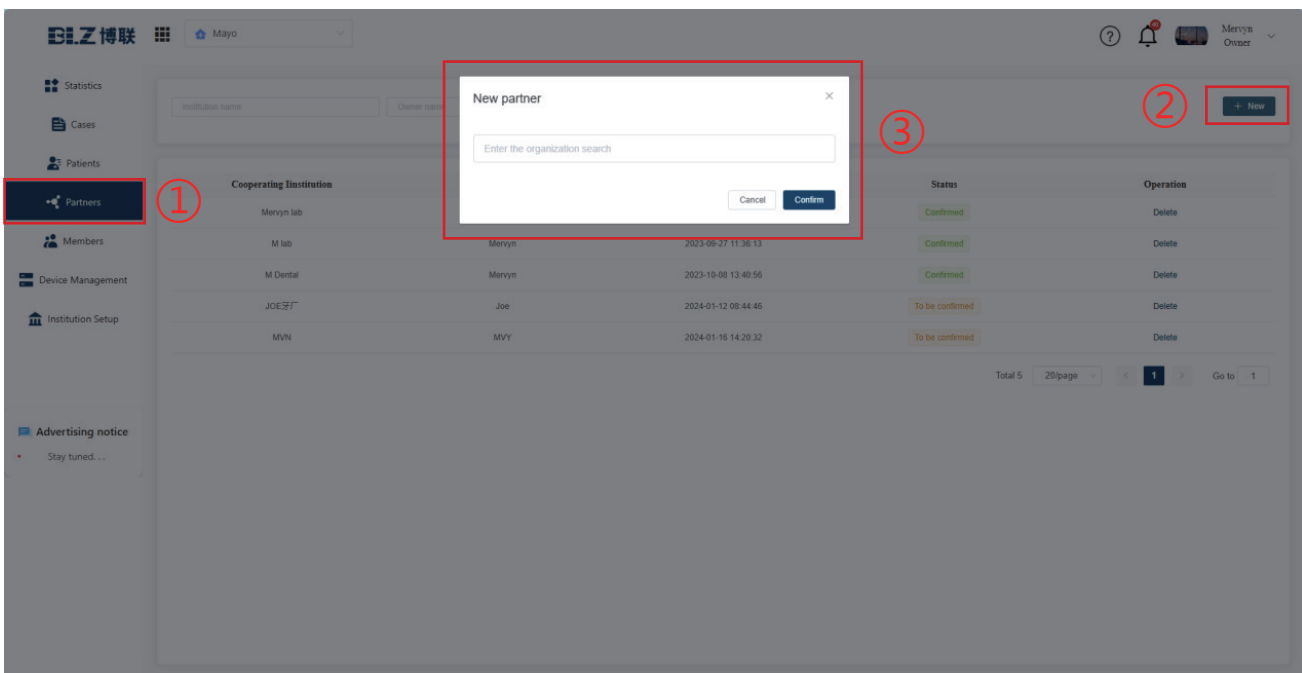


- ① Click “Partner”
- ② Click “New”
- ③ Input the name of dental lab and click “confirm”

- Dental lab logs in to agree the request.

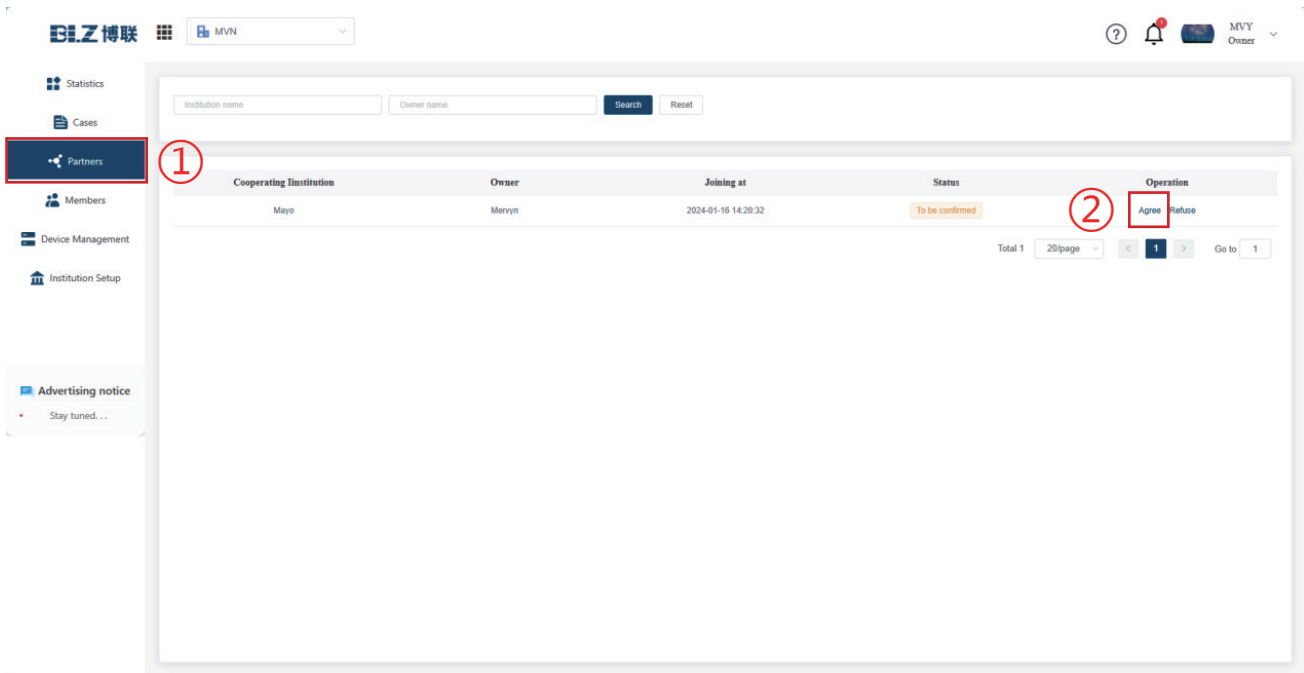

① After logging in, Click “Partners”

② Choose “agree”

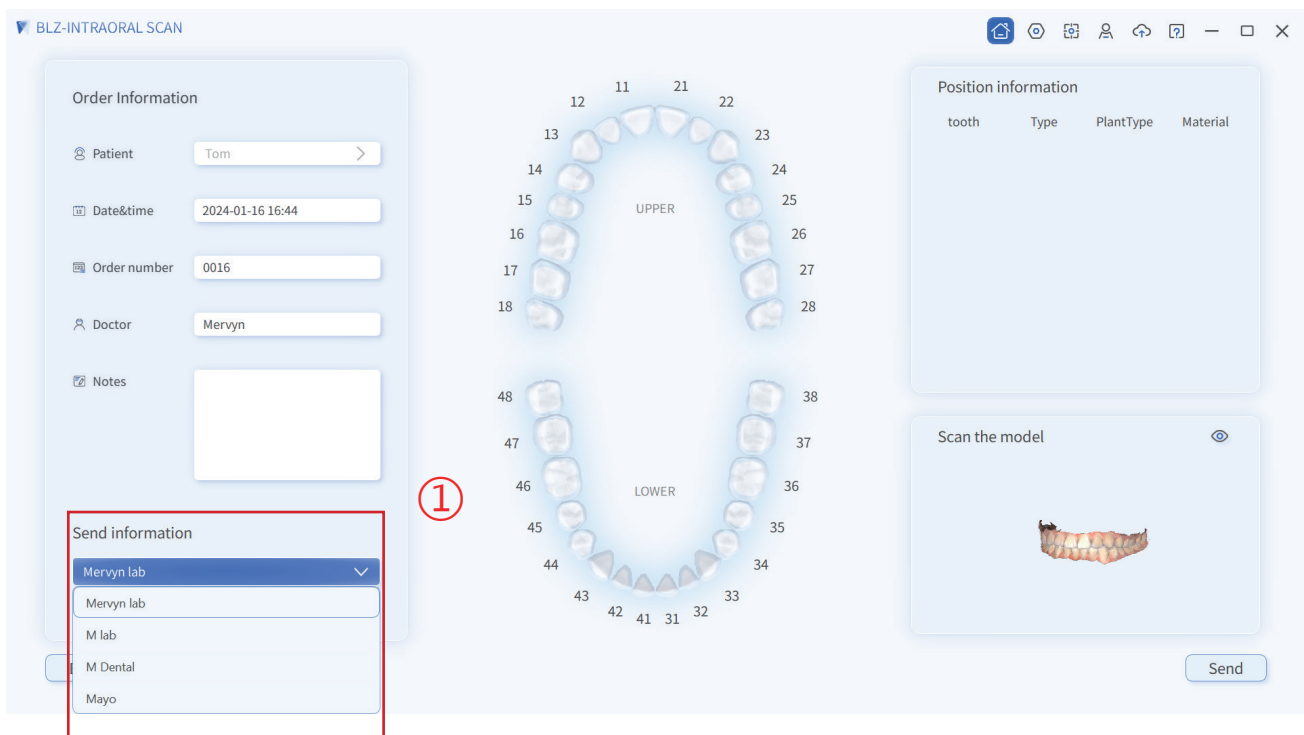


① After adding new partners, these dental labs can be selected when sending orders from scanning software.


3.Order management

3.1 Order management

- Upload orders from scanning software. Please make sure computer is connected to the internet.

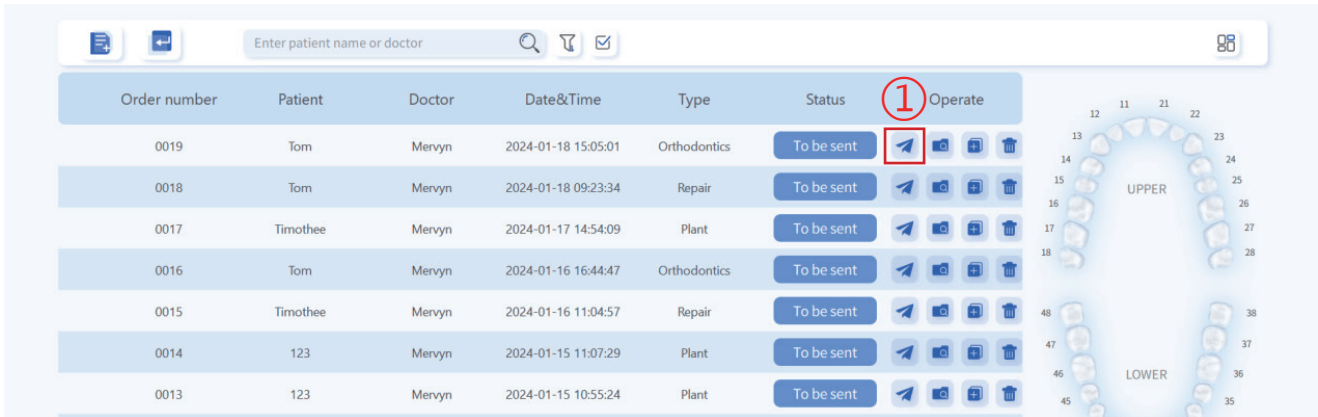


- ① Click to send orders to the BLZ cloud

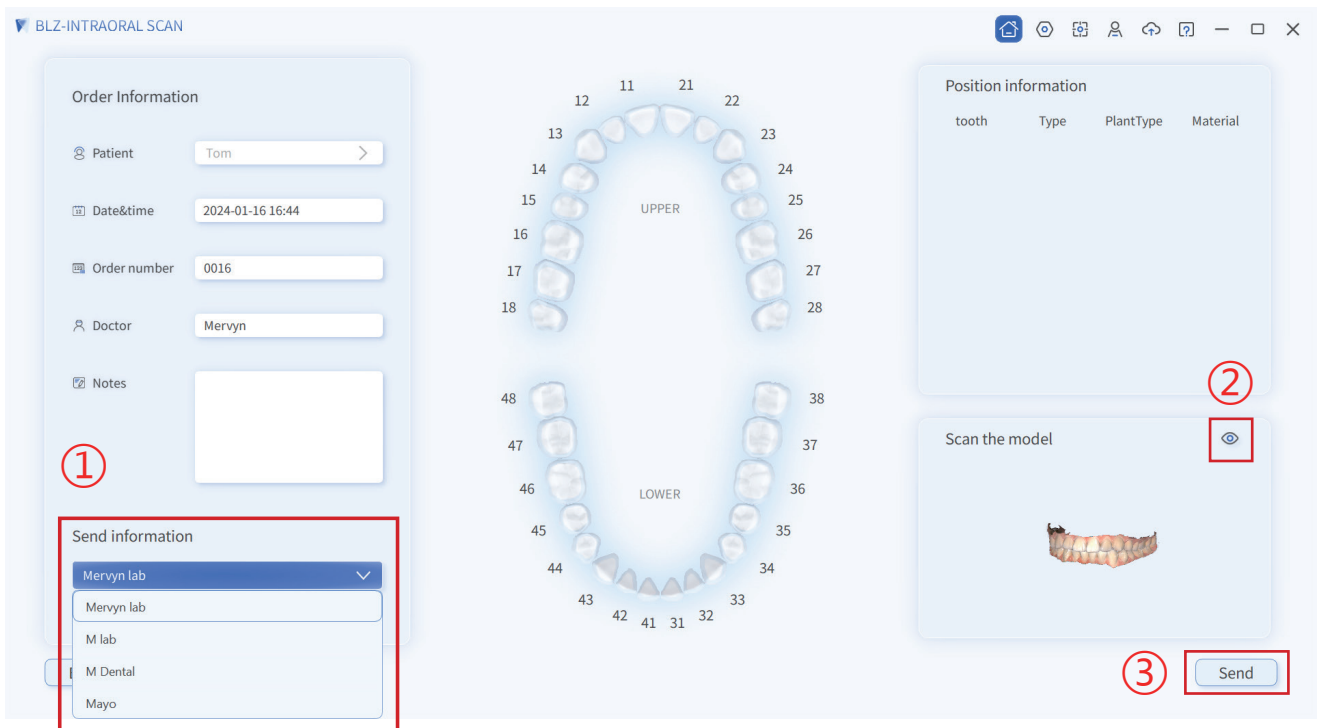

- ① Choose the dental lab that the order is sent to
- ② Click to preview the model
- ③ Click to send to BLZ cloud


3.2 Dental lab side

- After logging in, orders can be managed, including confirming orders, uploading design and uploading shipping info.

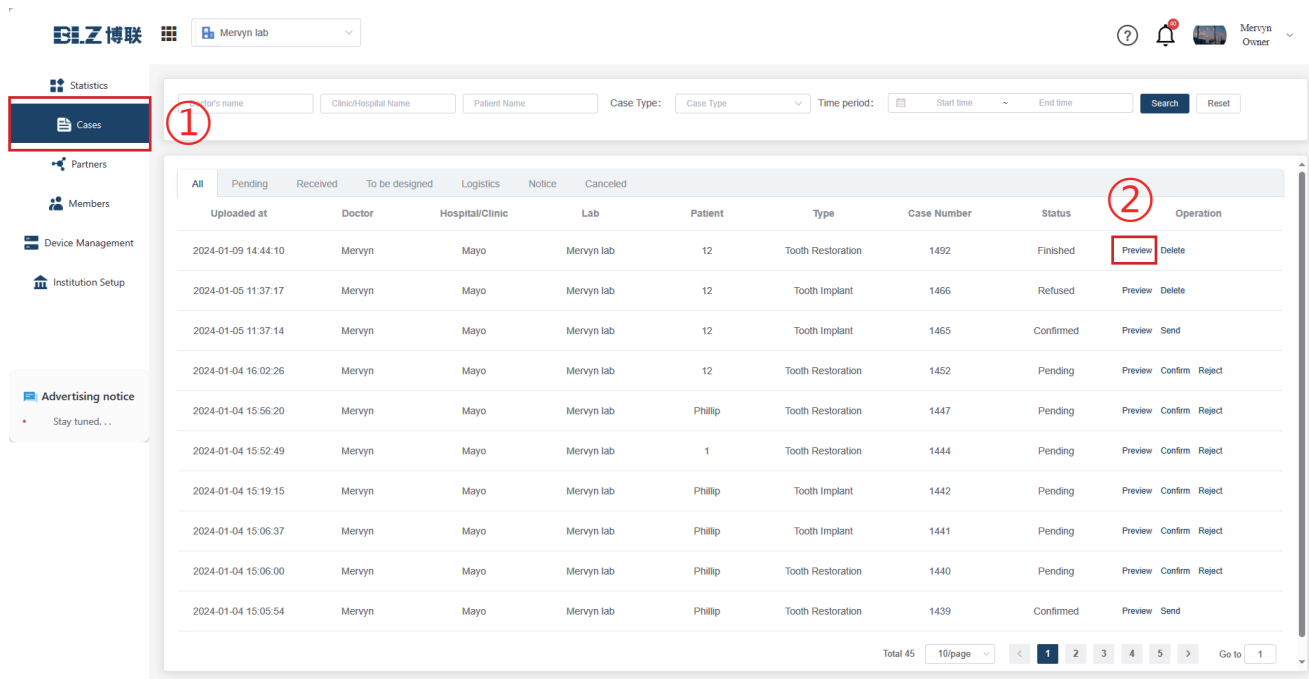

① After logging in, click “Cases”

② Click “Preview”

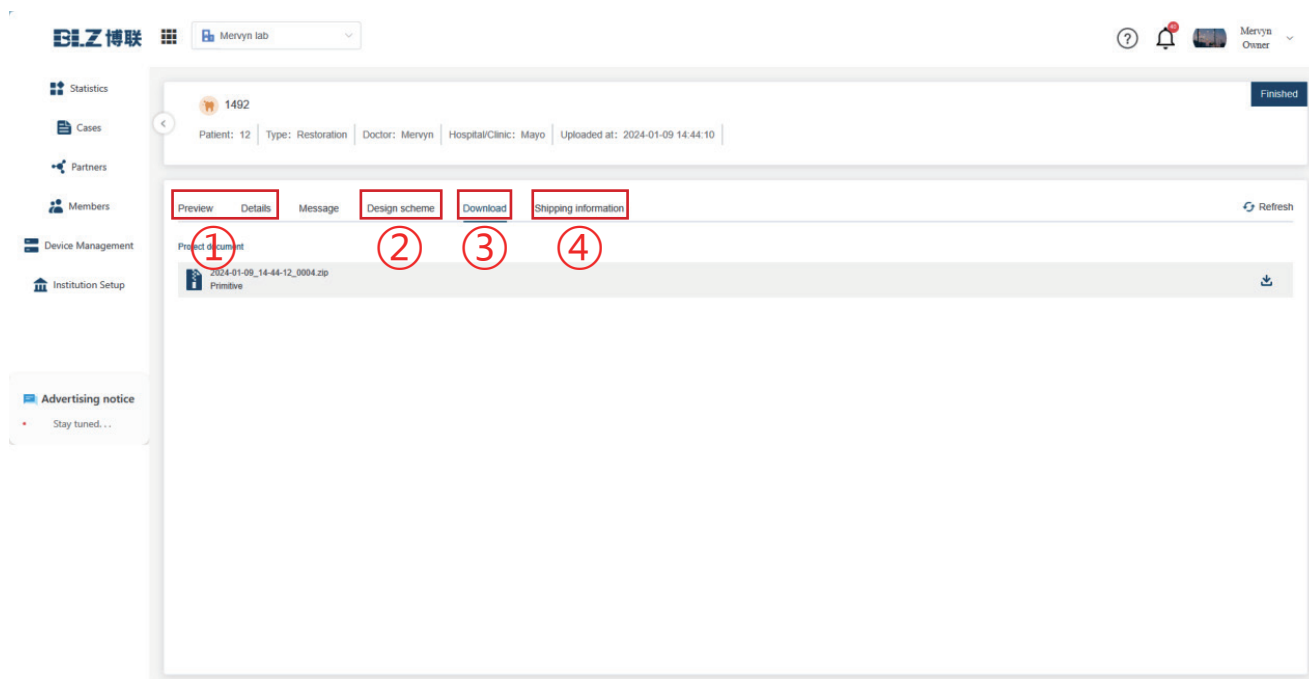

① Click to preview model and see details

② Click to upload design

③ Click to download model files

④ Click to upload shipping info



BLZ Technology (Wuhan) Co., Ltd.

Office Address:

C4166, Wuhan Overseas Scholar Business Park, No.11 Dongxin Road,
East Lake New Technology Development Zone, 430074, Wuhan, China


Manufactory Address:

3rd Floor, Building I, Tuochuang Industrial Park, No.15, Yangqiaohu
Avenue, Canglongdao Science and Technology Park, Jiangxia District,
430200, Wuhan, China


TEL: +86-27-6552 0255

WEB: www.blz-dental.com


Email: info@blz-dental.com

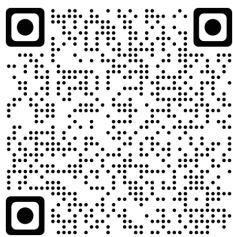

FACEBOOK Group

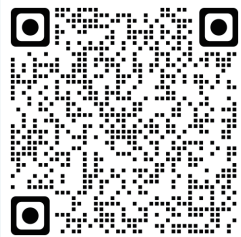

YOUTUBE

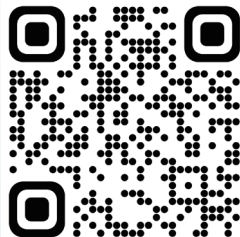

INS